
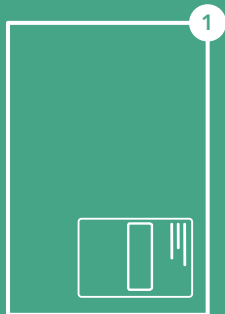


KIT COMPONENTS

1. Prepaid mailer
2. Specimen bag
3. Saliva kit
4. Shipping box
5. Barcode card

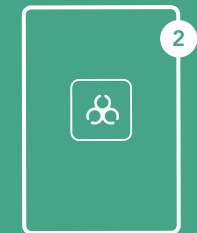

Do not remove plastic film from funnel cap.

The average time required to provide a good sample is approximately 2-5 minutes.

Be patient.

Before unscrewing the funnel from the tube, if your sample includes bubbles that reach the top of the funnel, let the sample settle to reduce potential spillage.

Seal the container immediately after providing your sample.

If you find it difficult to produce a complete sample, you may try to induce more saliva production by rubbing your cheeks or closing your mouth and wiggling your tongue. Placing a few grains of plain white sugar on your tongue will also make it easier to produce saliva.